

MiniBox+™

**Complete Pulmonary Function Testing.
Because Spirometry is just not enough.**




EASY

Easy to operate with minimum set-up. Fully automatic Lung Volumes and DLCO measurements. All results automatically generated.



ACCURATE

Clinically validated in over 500 patients at 8 medical centers, based on comparative data from leading body plethysmography devices.

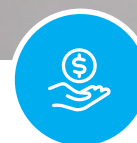

GREAT ROI

With all measurements fully reimbursable, you can achieve your ROI within months, with testing of only 2 to 5 patients daily.

Also available in MiniBox® configuration for spirometry and lung volumes only

MiniBox+

Complete, All-in-One Pulmonary Function Testing (PFT)

Desktop Plethysmography, Spirometry & Diffusion Testing

Easy to use, MiniBox+ provides accurate and repeatable, fully automatic measurements of Total Lung Capacity (TLC) and Residual Volume (RV), as well as Spirometry and Single Breath CO Diffusion (DLCO).

FULLY COMPLIANT WITH ATS/ERS GUIDELINES:

SPIROMETRY

Complete set of 30+ parameters including FEV1, FVC, SVC, MVV

LUNG VOLUMES

Same volumes as body plethysmography including TLC, FRC, RV

DIFFUSING CAPACITY


DLCO, VA, DLCO/VA
Single breath method

FULLY AUTOMATIC LVM AND DLCO

Manual measurement optional

ACCURATE & REPEATABLE

Automatically generated results

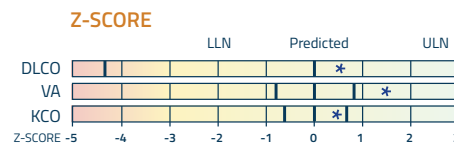

USER-FRIENDLY

Intuitive user interface, customizable reports, and a wide range of data export options

FULL SELECTION OF PREDICTED EQUATIONS

GLI reference values are standard for both Spirometry and Diffusion but easily customizable

PORTABLE

Lightweight with small footprint; Optional cart

PEDIATRIC USE

Measurements can be taken on children from age 5

GASLESS & CABINLESS PLETHYSMOGRAPHY

No sensation of claustrophobia; Can test obese, wheelchair-bound, or non-ambulatory patients

SELF-CALIBRATING

Due to unique sensor technology

CONNECTIVITY (OPTIONAL)

HL7 compliant, easily configurable, integrates with all major EHR systems

FeNO by  NIOX®

ADD-ON (OPTIONAL)

Fully integrated



The NIOX VERO device is attached to the MiniBox+ for transportation purposes only. FeNO tests should not be performed while attached to the MiniBox+.



* Reference values may not be available for DLCO for children under the age of 7

OTHER MEASUREMENTS

6MWT, Provocation

INNOVATIVE

Fully-automated Quality Control algorithms enable ease-of use and reliability

MiniBox+

Technical Specifications

General

Volume accuracy	±2%
Flow accuracy	±2.5% or 50 mL/s
Flow range	±14 L/s
Pressure range	±1 PSI, 7 kPa, Bi-directional
Dynamic resistance at 12L/s (spirometer)	<1.2 cmH ₂ O/(L/s)

Standards Compliance

Electrical	IEC 60601-1
EMC	IEC 60601-1-2

Physical

Full HD screen	10.8"
Dimensions	510x470x310 mm/ 20x18.5x12 in
Weight	8 kg / 18 lbs

Valve

Closing duration	100 ms
------------------	--------

Electrical

Power	100-240VAC, 50-60Hz
-------	---------------------

Gas Analyzer


Spectrometer	Infrared (CO, CH ₄)
Range	0-0.35%
Calibration	Automatic with quality control
Response time	t0-90 85 ms

With as little as 60 seconds of tidal breathing required for Lung Volume Measurements, this gasless, cabinless system provides maximum patient comfort.


EASE OF USE

- ✓ Automatic zeroing of analyzers before each measurement
- ✓ No stabilization time for plethysmography
- ✓ Warm-up time for analyzer less than 10 min
- ✓ Volume Gas graph for diffusion measurement
- ✓ Fast analyzer with response time t0-90 in less than 85 ms
- ✓ Time-to-test: as little as 20 minutes with pre- and post-spirometry

FLOW VOLUME

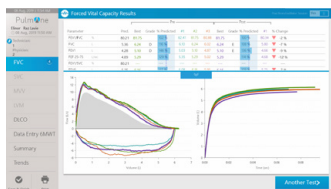

INCENTIVES

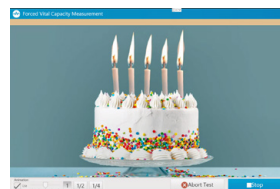


Complete mobility with dedicated MiniBox+ cart



For a complete list
of technical specs

MiniBox+ changes the economics of your practice

With MiniBox+, you can expand your medical practice while making Pulmonary Function Testing (PFT) more accessible to patients. Reimbursement for all testing measurements ensures a viable model for better patient care and better economics.

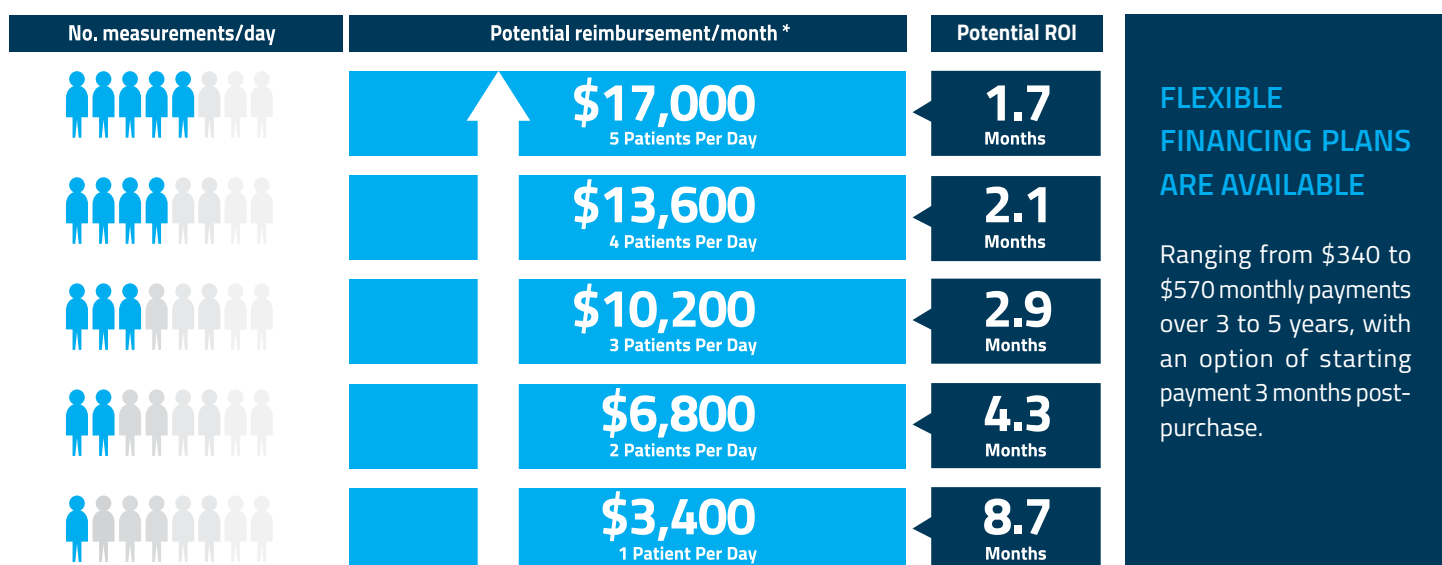

Consult with your financial advisor about the Section 179 tax benefits associated with the purchase of MiniBox+.

All tests are eligible for reimbursement*:

*subject to specific health plan policy

Code	Test
94010	Forced vital capacity - breathing capacity test
94060	Pre vs. Post bronchodilator - breathing capacity test
94200	Maximum breathing capacity
94726	Total lung capacity - determination of lung volume/plethysmography
94729	Diffusing Capacity
94375	Respiratory flow-volume loop
94618	Pulmonary stress testing (6-Minute-Walk-Test)
95012	Fractional exhaled nitric oxide (Add-on)

Coding should be based on medical necessity and is not intended to maximize reimbursement.

Expand your practice today. With MiniBox+.

ABOUT US

PulmOne is a global medical technology company dedicated to providing innovative solutions for complete Pulmonary Function Testing (PFT) in any clinical setting. At PulmOne, we develop patient-friendly, hassle-free, and budget-conscious devices that deliver a wide range of accurate measurements for the diagnosis, treatment, and monitoring of respiratory diseases.

PulmOne USA Tel: 1-844-PulmOne (1-844-785-6663)

Email us: usinfo@pulm-one.com

Visit us: www.pulm-one.com

PulmOne Advanced Medical Devices, Ltd. 14 Hamelacha St., Ra'anana 4366107 Israel

PulmOne
PFT Outside The Box